KIDYTED Kid's Growth Tracker

For safe use of the product, please read this User Manual first and save it for future reference. Please also keep safe the enclosed Warranty Card. Do not lose or discard it.

1 Product Use and Features

First, we would like to thank you for choosing Kidywolf's product. This product offers a more accurate and more reliable way to help you understand your kid's growth.

Product Application: weight measurement

Product Features:

1. Simple design – no button, no screen.
2. Auto on by stepping on.
3. Data displayed on and stored in the height ruler, no need to bend down to check.
4. Auto uploading of height ruler data after syncing with BonBaby APP, permanent Cloud storage.
5. Support kid's holding mode for kid's weight measurement.

2 Precautions


⚠ Precautions for Use

1. Do not jump on the product.
2. Leave the infrared transmitter uncovered while measuring.
3. Please place the scale on a flat and hard surface. Uneven surfaces or soft surfaces like carpet may result in inaccurate results.
4. Please use it in a room lower than 3.5m with a light-colored ceiling; otherwise, weight data may not show on the height ruler.
5. Please do not step on the scale with wet feet or when the scale surface is wet to avoid the risk of slipping.

⚠ Precautions for Storage

1. Do not soak it directly in water.
2. Do not leave it in high temperature, high humidity, dusty environment, direct sunlight or directly against the air conditioner.
3. Do not leave it in places where a strong impact may occur, vibration is present or on an uneven surface.


3 Product Illustration

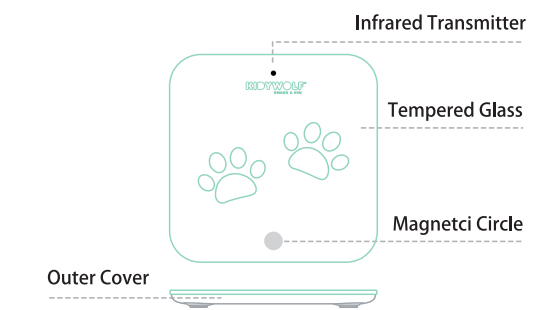The scale consists of tempered glass, an infrared transmitter, a magnetic circle, an outer cover, etc.

4 Accessories List

Please check that the following items are included in the package: a Kidyted Scale, 4 1.5V (AAA) batteries, and a User Manual. In case of anything missing, please contact our after-sales support.

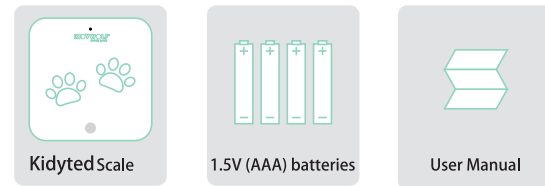5 Battery Installation Instructions

Battery installation instructions for Kidyted Scale:

1. Open the battery compartment cover.
2. Insert 4 AAA batteries into the battery compartment of the scale with the positive and negative ends facing the correct directions as indicated.
3. Put the battery cover back on in the reverse direction.


Friendly Reminders from Kidywolf:

1. If the Kidyted Scale is not going to be used for some time (3 months or above), please remove the batteries.
2. Do not mix old batteries with new ones or different types of batteries to avoid affecting the life of the product.
3. This product and used batteries shall be disposed in compliance with the electric product disposal regulations of the country of use.

6 How to Use


Step 1

Stand onto the scale with both feet and keep balance.

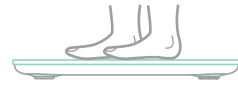Step 2

When the KIDYWOLF logo is flashing, the scale is measuring the weight.

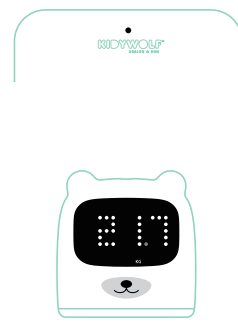Step 3

Once measurement is finished, the KIDYWOLF logo will remain on and weight data will appear on the height ruler.

7 Download BonBaby App

Search "BonBaby" on App Store or Google Play to download

Data will be synced with your APP through Bluetooth, enabling scientific and quantitative management of your baby's growth and offering professional analysis and advice.

Apple devices: iPhone 5 or above with version at least iOS 8.0, iPhone 4/4s excluded.
Android devices: 4.3 or above and supports Bluetooth 4.0/4.1

6 How to Use

Step 1

Turn and unfold the level ruler, move the ruler end closer to the scale and stick it onto the metal circle.

Step 2

Let your child stand upright on the scale and do not move. Pull the height ruler up, attach the level ruler onto the top of your child's head and keep it balanced.

Step 3

The measurement is completed after 3 seconds; the height data will be displayed on the LED dot matrix screen and flash 3 times.

7 Download BonBaby App

Search "BonBaby" on App Store or Google Play to download

Data will be synced with your APP through Bluetooth, enabling scientific and quantitative management of your baby's growth and offering professional analysis and advice.

Apple devices: iPhone 5 or above with version at least iOS 8.0, iPhone 4/4s excluded.
Android devices: 4.3 or above and supports Bluetooth 4.0/4.1

5 Battery Installation Instructions

Battery installation instructions for Kidyted Ruler:


1. Open the battery compartment cover.
2. Insert 2 AAA batteries into the battery compartment of the scale with the positive and negative ends facing the correct directions as indicated.
3. Put the battery cover back on in the reverse direction.

Low battery warning: if the battery icon flashes when turning the device on, it means that the batteries are running low. Please replace the batteries. Batteries with insufficient power may affect the accuracy of the measurement.

Friendly Reminders from NESPART SPRL:

1. If the Kidyted Ruler is not going to be used for some time (3 months or above), please remove the batteries.
2. Do not mix old batteries with new ones or different types of batteries to avoid affecting the life of the product.
3. This product and used batteries shall be disposed in compliance with the electric product disposal regulations of the country of use.

3 Product Illustration

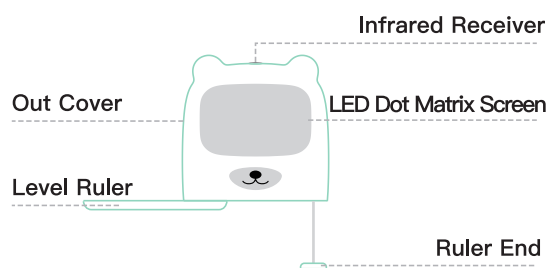The height ruler consists of an outer cover, an infrared receiver, a LED dot matrix screen, a level ruler and a ruler end.

2 Precautions

⚠ Precautions for Use

1. When measuring height, it is recommended to stick the ruler end to the middle of the metal circle on the scale.
2. Do not pull the measuring tape too fast so that the ruler end does not break away from the metal circle and hit the person on the scale.
3. When removing the ruler end from the magnetic circle, it is recommended to slide it out from the side of the metal circle.
4. After the measurement, please withdraw the measuring tape slowly. Withdrawing too fast may result in failure of placing it back to its original position.

⚠ Precautions for Storage

1. Do not soak it directly in water.
2. Do not leave it in high temperature, high humidity, dusty environment, direct sunlight or directly against the air conditioner.
3. Do not leave it in places where a strong impact may occur, vibration is present or on an uneven surface.

KIDYTED Kid's Growth Tracker

For safe use of the product, please read this User Manual first and save it for future reference. Please also keep safe the enclosed Warranty Card. Do not lose or discard it.

1 Product Use and Features

First, we would like to thank you for choosing Kidywolf's product. This product offers a more accurate and more reliable way to help you understand your kid's growth.

Product Application: height measurement

Product Features:

1. Original measuring solution for easy and accurate height measurement
2. Cute shape, small and portable, suitable for use at home
3. Integrated Bluetooth, complementary with BonBaby APP Auto transmission and permanent Cloud storage of data
4. Simple no-button design, turn and unfold the level ruler to activate height measuring
5. Off to sleep mode automatically, no need to turn off manually

8 Technical Specification

Product name:

KIDYTED Scale

Model:

KIDYTED-CY

Size:

312mm*312mm*23.5mm

Battery:

4 1.5V (AAA) batteries

Battery endurance:

6 months

Measurement:

3-150Kg

Division:

50g

Unit:

Kg/Lbs

Weight:

1.91Kg(without battery)

Operating conditions:

Temperature: 0°C~40°C

Humidity: ≤85%RH

Storage and transportation conditions:


Temperature: -20°C~60°C

Humidity: ≤85%RH

10 Maintenance and Care

1. Keep the Kidyted Scale away from heat.
2. This product is not waterproof. Do not let water drip or splash into the product. In case of exposure to water, please wipe with a dry cloth as soon as possible.
3. Wipe the scale surface with a clean and soft cloth. If necessary, wipe it with a wrung-out cloth and then dry it with a dry cloth.
4. Avoid wiping with a rough cloth in case of scratches on the scale.
5. Do not wipe this product with chemical or corrosive solvents.

11 Warranty and Service

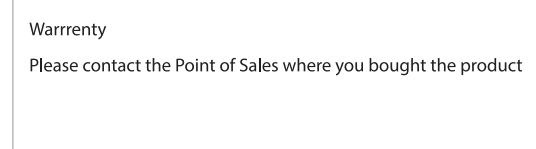Service:

1. 2-year warranty (from date of purchase) is included under proper usage.
2. Warranty doesn't cover parts other than Kidyted Scale itself.
3. Under these circumstances, repair will be charged:
 - A: Artificial destruction, inappropriate use, or used for commercial purposes
 - B: Warranty expired
 - C: Failure to present the Warranty Card
 - D: Content of the Warranty Card is inconsistent with the actual object, or blurred and unrecognizable
 - E: Dirt and scratches on the product itself or accessories.
 - F: Repaired and modified by the user, damaged by natural disaster or external force

12 Disclaimer and Copyright Statement

Disclaimer:

This manual is written in accordance with available information and contents are subject to change without notice. KIDYWOLF has made the best effort to guarantee the manual's accuracy and reliability. However, the company does not assume liability for loss or damage caused by omissions, inaccuracies or typographical errors in this manual. If there are discrepancies in this manual due to software upgrades, instructions and function dictated by the software.

Copyright Statement:

NESPART SPRL reserves all rights, and retains the final right to interpret and revise this handbook and the statements contained therein. Copy right of this manual belongs to NESPART SPRL. Without prior written permission, no company or individual is allowed to reproduce, copy or translate this manual or any information contained therein.

13 Multilingual Instruction download

- Download instruction manual in your language.
- Téléchargez le manuel d'instructions dans votre langue.
- Laden Sie die Bedienungsanleitung in Ihrer Sprache herunter.
- Scarica il manuale di istruzioni nella tua lingua.
- Descarga el manual de instrucciones en tu idioma.
- Download de handleiding in uw taal.

8 Technical Specification

Product name:

Kidyted Ruler

Model:

KIDYTED-CY

Size:

95mm*90mm*29mm

Battery:

2 1.5V (AAA) batteries

Measurement:

1-180cm

Wireless transmission:

BLE 4.1

Transmission distance:

10m

System requirements:

iOS: iPhone 5 or above with version at least iOS 8.0, iPhone 4/4s excluded.
Android: 4.3 or above and supports Bluetooth 4.0/4.1

Division:

0.2cm

Unit:

cm/in

Weight:

120g(without battery)

Battery endurance:

6 months

Operating conditions:

Temperature: 0°C~40°C

Humidity: ≤85%RH

Storage and transportation conditions:


Temperature: -20°C~60°C

Humidity: ≤85%RH

10 Maintenance and Care

1. Keep the Kidyted Ruler away from heat.
2. This product is not waterproof. Do not let water drip or splash into the product. In case of exposure to water, please wipe with a dry cloth as soon as possible.
3. Wipe the scale surface with a clean and soft cloth. If necessary, wipe it with a wrung-out cloth and then dry it with a dry cloth.
4. Avoid wiping with a rough cloth in case of scratches on the scale.
5. Do not wipe this product with chemical or corrosive solvents.

11 Warranty and Service

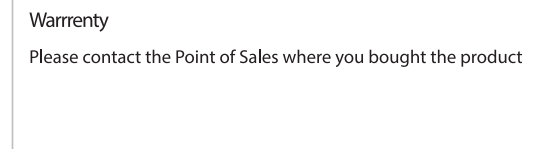Service:

1. 2-year warranty (from date of purchase) is included under proper usage.
2. Warranty doesn't cover parts other than Kidyted Scale itself.
3. Under these circumstances, repair will be charged:
 - A: Artificial destruction, inappropriate use, or used for commercial purposes
 - B: Warranty expired
 - C: Failure to present the Warranty Card
 - D: Content of the Warranty Card is inconsistent with the actual object, or blurred and unrecognizable
 - E: Dirt and scratches on the product itself or accessories.
 - F: Repaired and modified by the user, damaged by natural disaster or external force

12 Disclaimer and Copyright Statement

Disclaimer:

This manual is written in accordance with available information and contents are subject to change without notice. Kidywolf has made the best effort to guarantee the manual's accuracy and reliability. However, the company does not assume liability for loss or damage caused by omissions, inaccuracies or typographical errors in this manual. If there are discrepancies in this manual due to software upgrades, instructions and function dictated by the software.

Copyright Statement:

NESPART SPRL reserves all rights, and retains the final right to interpret and revise this handbook and the statements contained therein. Copy right of this manual belongs to NESPART SPRL. Without prior written permission, no company or individual is allowed to reproduce, copy or translate this manual or any information contained therein.


13 Multilingual Instruction download

- Download instruction manual in your language.
- Téléchargez le manuel d'instructions dans votre langue.
- Laden Sie die Bedienungsanleitung in Ihrer Sprache herunter.
- Scarica il manuale di istruzioni nella tua lingua.
- Descarga el manual de instrucciones en tu idioma.
- Download de handleiding in uw taal.

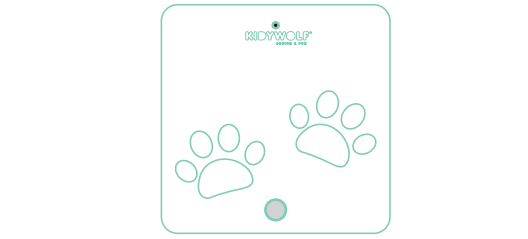Symptoms	Reasons	Solutions
LED screen doesn't light up when unfolding the level ruler.	Batteries are too low. Hardware is damaged.	Replace the batteries. Please contact us through www.kidywolf.com
Measurement data is not accurate.	Pull the height ruler too fast.	Please pull the height ruler slowly at a uniform speed.
Symbols	Instructions	
The screen shows the battery icon.		Low battery but will still work properly.
The screen shows the battery icon.		Batteries are too low to work properly.
The screen shows the icon "W".		Weight measurement mode.
The screen shows the icon "H".		Height measurement mode.
The screen shows the icon "OL".		Measuring weight exceeds the range.

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